

# Capp *Rondo*

## 4-PLACE PLATE SHAKER

**Capp *Rondo* 4-Place Plate Shaker** is the ideal choice for thorough, simultaneous mixing in up to 4 individual microplates of different types. The 3mm orbital motion movement guarantees the ultimate mixing performance and the brushless DC motor ensures long and maintenance free operation, making **Capp *Rondo* 4-Place Plate Shaker** the ideal device for extended runs.

**Capp *Rondo* 4-Place Plate Shaker** operates with the speed range from 200 to 1200 RPM. Its small and compact design requires minimum benchtop space, and the digital display gives the user an overview of the speed and time settings.


**Capp *Rondo* 4-Place Plate Shaker** can operate from 1 to 999 minutes, or without any time limits - in an infinite mode.

The rubber o-rings ensure an ultimate grip for all 4 plates. At the same time, loading and unloading is done effortlessly. Once loaded, all 4 plates can be easily accessed.

The pulse mode of **Capp *Rondo* 4-Place Plate Shaker** enables continuous mixing in clockwise and counterclockwise directions, changing in a time interval from 30 to 90 seconds. Thanks to the last run memory feature, the user can comfortably continue their work in the same mode as previously.

**Capp *Rondo* 4-Place Plate Shaker** is compatible with all commonly used skirted microplates and extremely easy to clean and maintain.



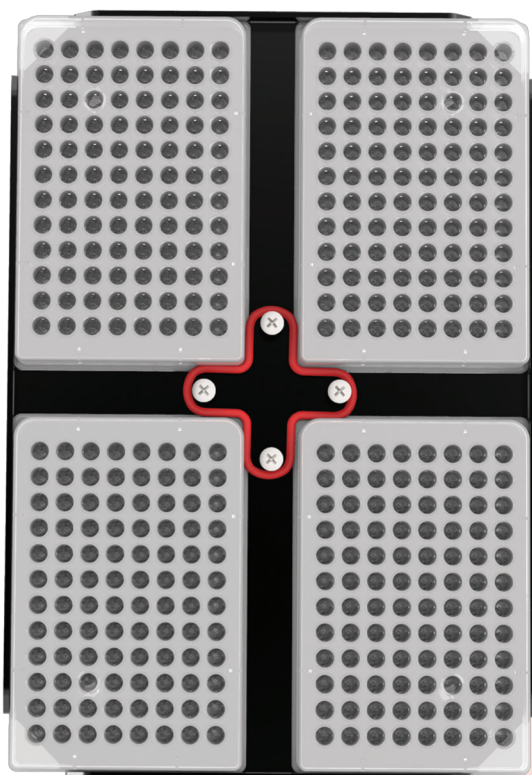

| Specifications:        | Capp <i>Rondo</i> 4-Place Plate Shaker |
|------------------------|----------------------------------------|
| Motor type             | BLDC motor                             |
| Orbital Motion         | 3mm                                    |
| Number of plates       | Up to 4                                |
| Variable Speed         | 200 rpm to 1200 rpm                    |
| Run Time               | 1 to 999 mins and infinity mode        |
| Pulse Mode             | Yes, Programmable (30 to 90 seconds)   |
| Weight                 | 5.4kg                                  |
| Maximum Volume         | 4 x 384 x 60µL                         |
| Inbuilt Counterbalance | Yes                                    |
| Dimensions (L x W x H) | 343 x 88 x 195mm                       |
| Digital Display        | Yes                                    |
| Input Voltage          | 24V                                    |
| Attachments            | Acc. All standard Skirted Microplates  |

# Ordering Information

| Cat No.  | Description                             |
|----------|-----------------------------------------|
| CRP-412X | CappRondo 4-Place Plate Shaker 1200 rpm |

